

Lab Photo Reactor Option

- Compatible with the Lab Reactor system
- Tunable LED irradiation source, with 6 different wavelengths
- LED lighting intensity higher than 100 mW/cm²
- Efficient light penetration with both sides of glass fluidic module illuminated
- Extended LED lifetime due to efficient liquid cooling
- Wireless control of wavelength selection and intensity



CORNING

The future flows through Corning® Advanced-Flow™ Reactors

Lab Reactor

CORNING

EMEA and NSA

Corning S.A.S.
Reactor Technologies
7 bis Avenue de Valvins
CS 70156 Samois sur Seine
77215 Avon Cedex, FRANCE
ph. +33 1 64 69 71 07
fax +33 1 64 69 70 59
reactors@corning.com

CHINA

Corning China (Shanghai) RHQ
No. 358 Lu Qiao Road
Jinqiao Export Processing Zone, Pudong
Shanghai 201206, CHINA
ph. +86 21 22152888 *1408
fax +86 21 621522988
reactor.asia@corning.com



3265 Sunset Lane, Hatboro, PA 19040
www.sentinelprocess.com
P: 800.345.FLOW F: 888.329.9669

www.corning.com/reactors


© 2018 Corning Incorporated. All rights reserved

Lab Reactor

Stepping into flow chemistry

Features

- Plug and play open system including data monitoring
- Ready to start, with thermostat integrated
- Up to 3 liquid lines, 1 gas line and a back pressure regulator for pressure control
- High chemical durability due to a full metal free system
- Outstanding mixing and heat exchange with patented HEART design
- Low internal volume
- Seamless scale-up with other Advanced-Flow™ Reactor products

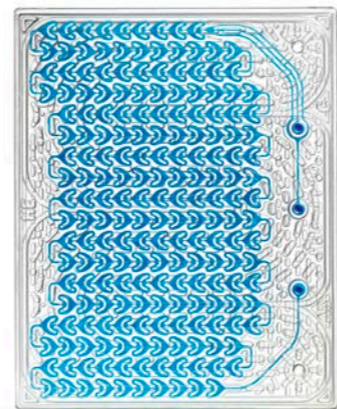

Fluidic module size:
155 x 125 mm



Size:
45 x 48 x 52 cm
(L x W x H)

Size:
40 x 38 x 45 cm
(L x W x H)

Technical Specifications

FLOW RATE	TEMPERATURE	PRESSURE	MATERIALS	FLUIDIC MODULE	OPTIONS
2 to 10 ml/min	-40°C to 200°C	Up to 18 barg	Glass PFA / PTFE Perfluoroelastomer	2.7 ml internal volume	3 rd liquid dosing line 2 nd glass fluidic module Lab Photo Reactor module

Mass Transfer 100 x better *

Heat Transfer 1000 x better *

Reaction Volume 1000 x lower *

Residence Time Distribution 50 x better *

* compared to batch reactors